




# Rocket<sup>®</sup> Infertility Range



**Rocket**medical

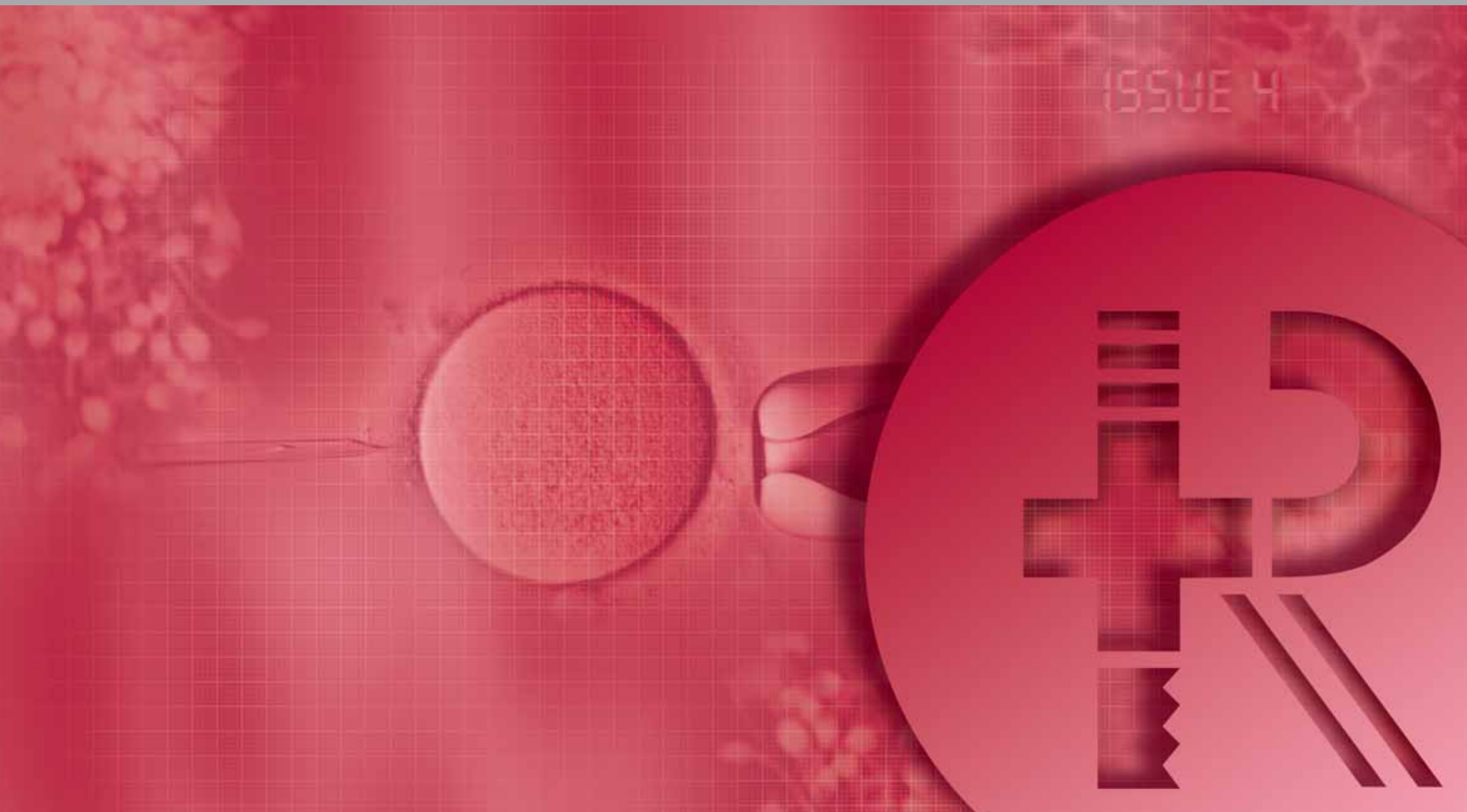

# Rocket® HSG & HyCoSy™ Catheter



The **Rocket® HSG & HyCoSy™** Catheter range provides a reliable, comfortable and cost effective method of performing uterine canalisation for hysterosalpingography and hysterosalpingo-contrast sonography.

## NEW INTEGRAL WIRE STYLET

- Gives improved malleability without using tiresome removable stylets.
- Improves ease of insertion.
- Increased value for money.
- CONTRAST PATH TAP ensures retention of contrast when further volumes of media are required.
- STERILE - Packed in peel pouch with syringe attached, cartons of 10 catheters. Intended for one time use.


Our newest HSG & HyCoSy™ Catheters also have all the trusted features of our regular catheter:

- LATEX FREE balloon minimises risk of patient tissue reactions.
- SMOOTH TIP & LOW PROFILE BALLOON ensures low trauma insertion and removal - ideal for the nulliparous cervix.
- REMOVABLE STIFFENING SHEATH to aid insertion in the heavily retroverted or antverted cervix.
- COMPATIBLE with all current contrast media.

Ordering Information	Pack	Code
<b>Rocket® HSG &amp; HyCoSy™ Catheter</b> Length 280mm, pre-set syringe supplied 1.5ml Maximum balloon volume: 3ml	5FG	R57422
	7FG	R57424
<b>Rocket® HSG &amp; HyCoSy™ Catheter with integral wire stylet</b> Length 280mm, pre-set syringe supplied 1.5ml Maximum balloon volume: 3ml	5FG	R57422-00-ST
	7FG	R57424-00-ST

### Specifications:

OD: 5FG & 7FG  
 Length: 280mm active length.  
 Balloon: Latex free, supplied with pre-set syringe: 1.5ml  
 Maximum balloon volume: 3ml  
 Sterilisation: Gamma irradiated @ 25kGy/2.5 MRads.

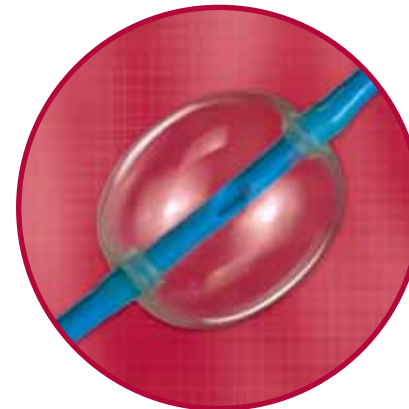



**Rocket**medical



## rocketmedical.com

### **UK Customer Services**

Rocket Medical plc  
Sedling Road, Washington  
Tyne and Wear  
NE38 9BZ England  
Telephone: +44 (0) 191 419 6988  
Fax: +44 (0) 191 419 6989  
Email: [homesales@rocketmedical.com](mailto:homesales@rocketmedical.com)

### **Head Office & Export Customer Services**

Rocket Medical plc  
Imperial Way  
Watford, Herts  
WD24 4XX England  
Telephone: +44 (0) 1923 651404  
Fax: +44 (0) 1923 240334  
Email: [exportsales@rocketmedical.com](mailto:exportsales@rocketmedical.com)

### **US Office**

Rocket Medical  
150 Recreation Park Drive, Unit 3  
Hingham, MA 02043, USA  
Customer Services: 1-800-707-ROCKET (7625)  
Telephone: 1-781-749-6223  
Fax: 1-781-749-6235  
Email: [usa@rocketmedical.com](mailto:usa@rocketmedical.com)